SONY. DIGITAL GRAPHIC PRINTER UP-D897

SERVICE MANUAL 1st Edition

∧警告

このマニュアルは、サービス専用です。

お客様が,このマニュアルに記載された設置や保守,点検,修理などを行うと感電や火災, 人身事故につながることがあります。

危険をさけるため、サービストレーニングを受けた技術者のみご使用ください。

⚠ WARNING

This manual is intended for qualified service personnel only.

To reduce the risk of electric shock, fire or injury, do not perform any servicing other than that contained in the operating instructions unless you are qualified to do so. Refer all servicing to qualified service personnel.

⚠ WARNUNG

Die Anleitung ist nur für qualifiziertes Fachpersonal bestimmt.

Alle Wartungsarbeiten dürfen nur von qualifiziertem Fachpersonal ausgeführt werden. Um die Gefahr eines elektrischen Schlages, Feuergefahr und Verletzungen zu vermeiden, sind bei Wartungsarbeiten strikt die Angaben in der Anleitung zu befolgen. Andere als die angegeben Wartungsarbeiten dürfen nur von Personen ausgeführt werden, die eine spezielle Befähigung dazu besitzen.

AVERTISSEMENT

Ce manual est destiné uniquement aux personnes compétentes en charge de l'entretien. Afin de réduire les risques de décharge électrique, d'incendie ou de blessure n'effectuer que les réparations indiquées dans le mode d'emploi à moins d'être qualifié pour en effectuer d'autres. Pour toute réparation faire appel à une personne compétente uniquement.

Table of Contents

Manual Structure	3-1-1. Outline
	3-1-2. Block Diagram of Electrical Block3-1
Purpose of this manual3	3-2. Head/Mechanical Control Block3-1
Related manuals3	3-2-1. Outline
	3-2-2. Platen Motor Control3-1
	3-2-3. Thermal Head UP/DOWN Control3-1
1. Service Overview	3-2-4. Monitoring of Door Sensor3-2
	3-2-5. Monitoring of Paper Sensor3-2
1-1. Board Location1-1	3-2-6. Monitoring of Head Temperature Sensor 3-2
1-2. Main Parts Location1-1	3-2-7. Control of Head Fan Motor (for Head Cooling) 3-2
1-3. Removing/Installing the Cabinet	3-2-8. Read of Buttons3-3
1-3-1. Top Cover1-1	3-3. USB Interface
1-3-2. Front Panel Subassembly1-2	3-4. Thermal Head Block
1-3-3. Rear Panel1-2	3-4-1. Structure
1-4. Replacing the Main Parts1-3	3-4-2. Basic Operation
1-4-1. Switching Regulator/DC Fan1-3	3-4-3. Temperature Correction
1-4-2. Mechanical Deck Assembly1-3	3-4-4. Correction of Resistor Count3-4
1-4-3. Thermal Head1-4	
1-4-4. Stepping Motor/DC Motor1-5	
1-4-5. Cam Shaft Assembly1-6	4. Troubleshooting
1-4-6. Pinch Arm Assembly1-7	
1-4-7. Platen Roller1-7	4-1. Print result is not satisfactory4-1
1-4-8. Timing Belt1-8	4-2. Print result density is too high or low4-1
1-5. Procedure Required for Replacement1-9	4-3. Trouble of determining presence or absence of paper 4-2
1-6. Unleaded Solder1-9	4-4. Thermal head UP/DOWN operation trouble4-2
	4-5. Feed operation trouble4-3
	4-6. Trouble of determining door open/close4-3
2. Electrical Alignment	
2-1. Procedure before and after Adjustment2-1	5. Service Mode (Self-diagnosis Function)
2-1-1. Procedure of Printing the User Setting Record	
(Before Adjustment)2-1	5-1. Startup Procedure5-1
2-1-2. User Setting Check Procedure	5-2. Service Mode Menu5-1
(After Adjustment)2-1	5-3. Test Pattern Printing5-1
2-1-3. Factory Setting2-1	5-4. Test Pattern Print Count5-3
2-2. Calibration and Electrical Conductivity Check2-2	5-5. LCD Lighting Check5-3
2-3. Head Voltage Adjustment2-3	5-6. LED Lighting Check
2-4. Initialization of Print Count History2-4	5-7. Front Panel Volume Calibration5-4
2-5. Firmware Version Upgrade2-5	5-8. Log Printing
3. Circuit Description	5-10. Temperature Indication
	5-11. FEED Operation
3-1. System Control Block	5-12. Menu Lock

6. Spare Part	S
---------------	---

6-1.	Notes on Repair Parts			
6-2.	Exploded Views	***********	**********	6-2
6-3.	Electrical Parts List			6-7
6-4.	Supplied Accessories			
7.	Block Diagram			
	Overall			7-1
_	Oshamatia Diagram			
8.	Schematic Diagram	IS		
	KY-572			8-7
	MA-132			
	SE-768			
	SE-769			
	SU-112			
	Frame Wiring			
	Frame Wiring	*************		0-10
9.	Board Layouts			
-	,	•		
	KY-572			9-1
	SE-768			
	SE-769			
	SU-112			
				9-2

Manual Structure

Purpose of this manual

This manual is the service manual of Digital Graphic Printer UP-D897.

This manual describes the information on maintenance and the service information such as service overview, electrical alignment, circuit description, troubleshooting, and service mode (self-diagnosis function).

Related manuals

In addition to this "Service Manual", this unit is provided with the manual below.

 "Operating Instruction" PDF (Included in the CD-ROM Supplied for products.)

Part No.: 3-863-306-0X


These manuals describes the information required for the actual management and operation of this unit.

"Semiconductor Pin Assignments" CD-ROM (Available on request)
 This "Semiconductor Pin Assignments" CD-ROM allows you to search for semiconductors used in B&P Company equipment.
 Part number: 9-968-546-XX


3

Section 1 Service Overview

1-1. Board Location

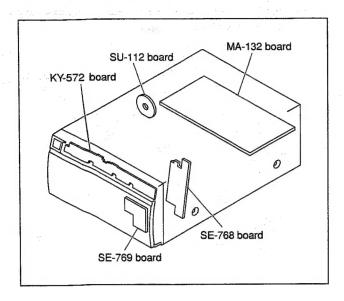


1-2. Main Parts Location

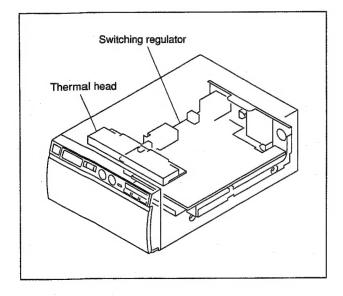

1-3. Removing/Installing the Cabinet

1-3-1. Top Cover


1. Remove the five M3 case securing screws, then remove the top cover.

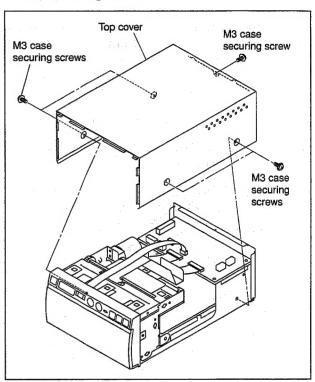


2. Attach the top cover in the reverse order.

1-3-2. Front Panel Subassembly

- 1. Remove the top cover. (Refer to Section 1-3-1.)
- Release the lock arm of the mechanical deck block (bottom surface) in the direction of the arrow, then open the door panel.
- 3. Remove the four hooks in both directions of the arrows, then remove the door panel in the direction of the arrow.
- 4. Disconnect the flexible flat cable from the connector (CN201) on the MA-132 board.

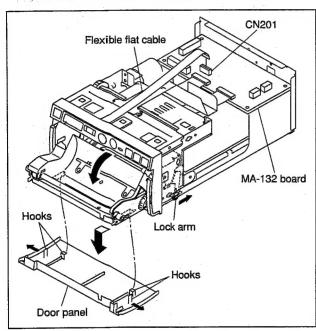


5. Remove the two screws and five hooks, then remove the front panel subassembly.

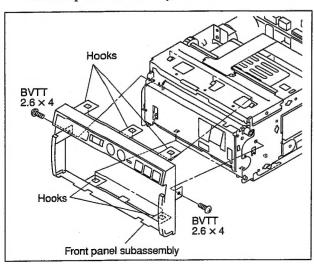

 Attach the front panel subassembly in the reverse order of steps 1 to 5.

1-3-3. Rear Panel

- 1. Remove the top cover. (Refer to Section 1-3-1.)
- 2. Remove the two screws (PS3 \times 6) and the screw (B3 \times 6), then remove the rear panel.

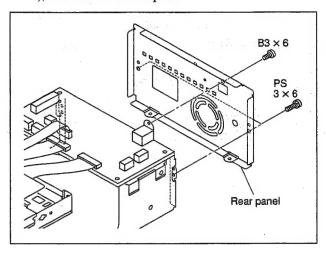

3. Attach the rear panel in the reverse order of steps 1 and 2.


1-4. Replacing the Main Parts

1-4-1. Switching Regulator/DC Fan

Note

There are two types of DC fan mounting connectors. Connect the harness of the DC fan for power supply to the connector (CN601). The connector (CN603) is provided for mounting the optional DC fan for cooling the thermal head.


- 1. Remove the top cover. (Refer to Section 1-3-1.)
- 2. Remove the front panel subassembly. (Refer to Section 1-3-2.)
- 3. Remove the rear panel assembly. (Refer to Section 1-3-3.)
- 4. Remove the mechanical deck assembly. (Refer to Section 1-4-2.)
- 5. Disconnect all harnesses and flexible flat cables from the MA-132 board.
- 6. Remove the four screws (BVTT2.6 × 6), then remove the MA-132 board.
- 7. Remove the two screws (PS3 \times 6), then remove the shield plate and insulating sheet (MA).
- 8. Remove the two screws (P3 \times 10), then remove the DC fan.
- 9. Remove the four screws (PS3 × 6), then remove the switching regulator.

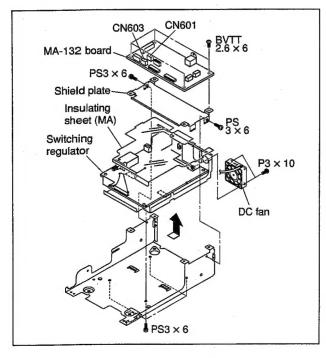

10. Attach the switching regulator in the reverse order of steps 1 to 9.

1-4-2. Mechanical Deck Assembly

- 1. Remove the top cover. (Refer to Section 1-3-1.)
- 2. Remove the front panel subassembly. (Refer to Section 1-3-2.)
- Remove the switch rod retainer, then remove the power switch rod.
- 4. Disconnect the three harnesses from the connectors (CN204, CN205 and CN602) on the MA-132 board.
- Disconnect the two flexible flat cables from the connectors (CN206 and CN302) on the MA-132 board.
- 6. Remove the three screws, then remove the mechanical deck assembly.

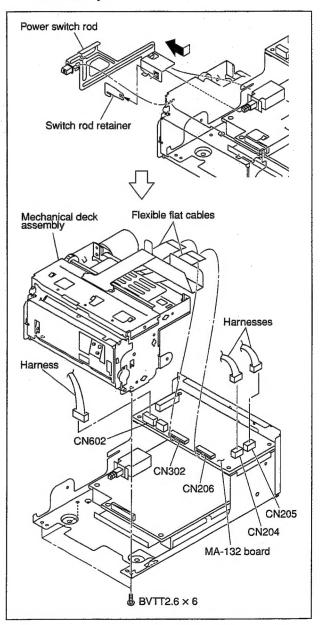

7. Attach the mechanical deck assembly in the reverse order of steps 1 to 6.

1-4-3. Thermal Head

Note


After replacing the thermal head, perform the head voltage adjustment. For the adjustment procedure, refer to Section 2-3.

- 1. Remove the top cover. (Refer to Section 1-3-1.)
- 2. Remove the front panel assembly. (Refer to Section 1-3-2.)
- 3. Remove the mechanical deck assembly. (Refer to Section 1-4-2.)
- Remove the two screws (BVTT2.6 × 6), then remove the MD top chassis and three head springs.

 Note

When removing the MD top chassis, be sure to remove it slowly in the direction of the arrow A to prevent the head springs from popping out.

- 5. Release the lock arm in the direction of the arrow, then open the MD door chassis.
- 6. Remove the heat sink in the direction of the arrow B.
- 7. Remove the two screws (PSW3 × 8), then remove the thermal head.
- 8. Disconnect the two flexible flat cables from the thermal head.

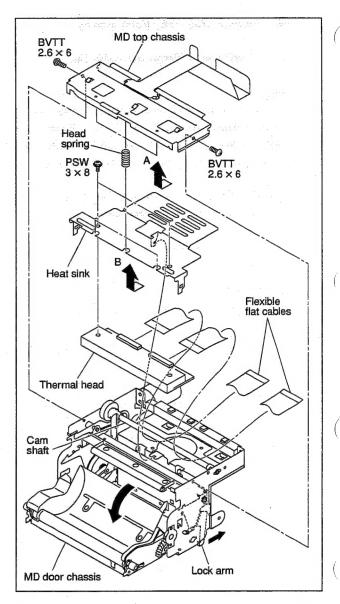

9. Attach the thermal head in the reverse order of steps 1 to 8.


Note

When connecting the flexible flat cables, route them under the cam shaft.

1-4-4. Stepping Motor/DC Motor


Removal

- 1. Remove the top cover. (Refer to Section 1-3-1.)
- 2. Remove the front panel assembly. (Refer to Section 1-3-2.)
- 3. Remove the mechanical deck assembly. (Refer to Section 1-4-2.)
- 4. Remove the helical extension spring from the tension arm assembly.
- 5. Remove the four screws (BVTT2.6 × 6), then remove the stepping motor bracket.
- 6. Remove the two screws (PS3 × 4), then remove the DC motor.
- 7. Remove the two screws (PS3 \times 6), then remove the stepping motor.

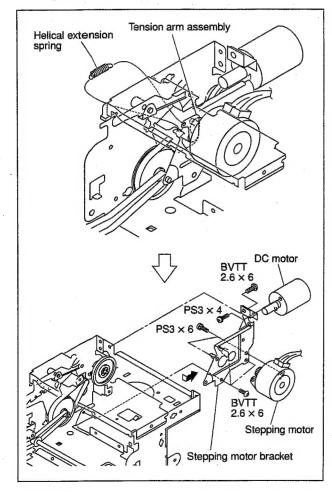

Installation

- 8. Attach the DC motor with the two screws.
- Insert the head up lever shaft (L) into the hole [A] and insert the speed reducing gear shaft (L) into the hole [B], then attach with the three screws.
- 10. Hang the timing belt on the stepping motor gear.

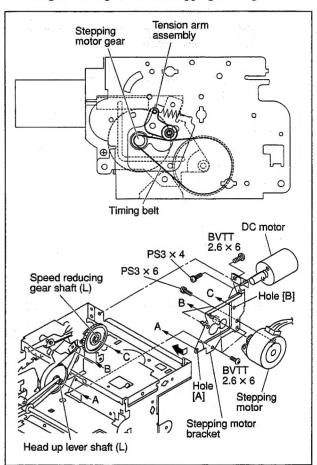


11. Assemble this unit in the reverse order of steps 1 to 4.

1-4-5. Cam Shaft Assembly


Removal

- 1. Remove the top cover. (Refer to Section 1-3-1.)
- 2. Remove the front panel assembly. (Refer to Section 1-3-2.)
- 3. Remove the mechanical deck assembly. (Refer to Section 1-4-2.)
- 4. Remove the thermal head. (Refer to Section 1-4-3.)
- 5. Remove the paper holder assembly. (Refer to step 4 of Section 1-4-7.)
- 6. Remove the two bearings by rotating them in the direction of the arrows A.
- 7. Remove the cam shaft assembly in the direction of the arrow B.
- 8. Remove the cam (R), cam gear, cam (L), noise reducing spring (30) and three parallel pins from the cam shaft.

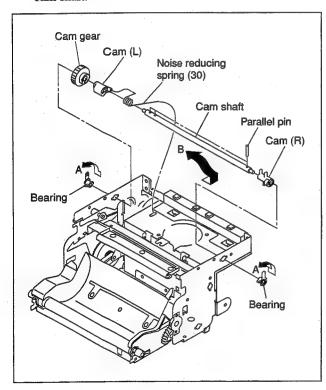

Installation

- 9. Attach the cam shaft assembly in the reverse order of steps 6 to 8.
- 10. Insert the portion A of noise reducing spring (30) into the hole of chassis, then insert the portion B into the hole of cam (L) by rotating it in the direction of the arrow.

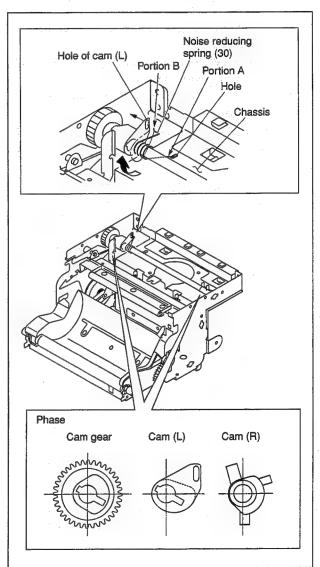


11. Assemble this unit in the reverse order of steps 1 to 5.

1-4-6. Pinch Arm Assembly

Note

When replacing the pinch arm subassembly, be extremely careful not to touch it with bare hands or dirty gloves. If the pinch arm subassembly is contaminated, clean it with the ethyl alcohol.

- 1. Remove the top cover. (Refer to Section 1-3-1.)
- 2. Remove the front panel assembly. (Refer to Section 1-3-2.)
- Remove the mechanical deck assembly. (Refer to Section 1-4-2.)
- 4. Remove the thermal head. (Refer to Section 1-4-3.)
- Remove the paper holder assembly.
 (Refer to step 4 of Section 1-4-7.)
- 6. Remove the three helical extension springs.
- Remove the four bearings, then remove the door lock assembly and pinch arm assembly in the direction of the arrow.
- 8. Remove the pinch arm assembly from the door lock assembly.

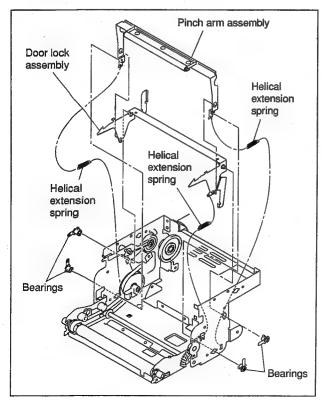

9. Attach the pinch arm assembly in the reverse order of steps 1 to 8.


After attaching the pinch arm assembly, clean it with ethyl alcohol.

1-4-7. Platen Roller


Note

When replacing the platen roller, be extremely careful not to touch it with bare hands or dirty gloves. If the roller is contaminated, clean it with the ethyl alcohol.

- 1. Remove the top cover. (Refer to Section 1-3-1.)
- 2. Remove the front panel assembly. (Refer to Section 1-3-2.)
- 3. Remove the mechanical deck assembly. (Refer to Section 1-4-2.)
- Remove the two screws, then remove the paper holder assembly.
- 5. Remove the helical extension spring from the tension arm assembly.

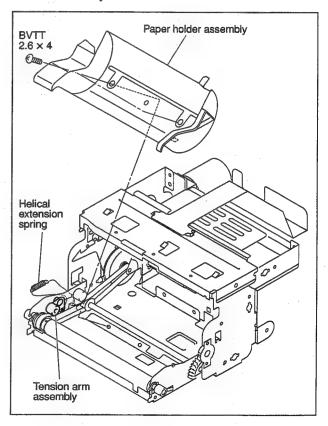


- 6. Remove the two platen bearings by rotating them in the direction of the arrow A.
- 7. Remove the platen roller in the direction of the arrow B.
- Remove the platen pulley and parallel pin from the platen roller.

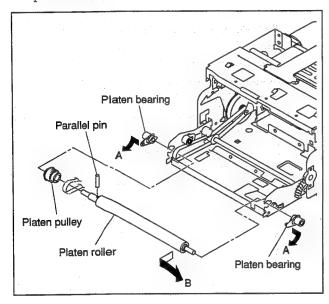

9. Attach the platen roller in the reverse order of steps 1 to 8.

1-4-8. Timing Belt

- 1. Remove the top cover. (Refer to Section 1-3-1.)
- 2. Remove the front panel assembly. (Refer to Section 1-3-2.)
- 3. Remove the mechanical deck assembly. (Refer to Section 1-4-2.)
- 4. Remove the paper holder assembly. (Refer to step 4 of Section 1-4-7.)
- 5. Remove the two screws, then remove the oil damper.
- 6. Remove the damper gear.
- Remove the door releasing spring from the hole of MD door chassis, then remove the MD door chassis assembly.
- 8. Remove the fulcrum spacer, then remove the door releasing spring.
- 9. Pull out the fulcrum shaft from the MD door chassis assembly in the direction of the arrow.
- 10. Remove the timing belt [A] and timing belt [B].

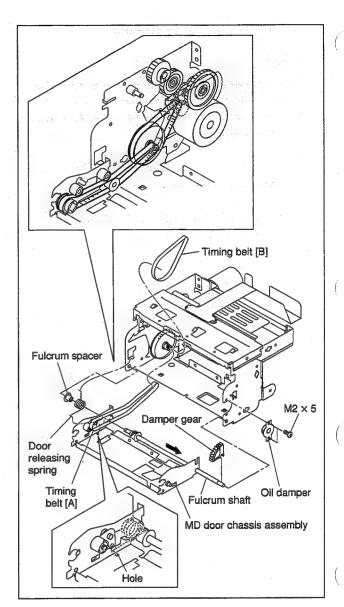

11. Attach the timing belt [A] and timing belt [B] in the reverse order of steps 1 to 10.

1-5. Procedure Required for Replacement

When replacing the following parts, perform each required process according to the table below. For the procedure before and after performing the required process, refer to Section 2-1.

Parts	Required process	Reference
MA-132 board	Firmware Version Upgrade	Section 2-5
IC103 (MA-132 board)		
KY-572 board	Calibration and Electrical Conductivity Check	Section 2-2
Switching regulator	Head Voltage Adjustment	Section 2-3
Thermal head		

1-6. Unleaded Solder

Boards requiring use of unleaded solder are printed with a lead free mark (LF) indicating the solder contains no lead.

(Caution: Some printed circuit boards may not come printed with the lead free mark due to their particular size.)



Notes

- Be sure to use the unleaded solder for the printed circuit board printed with the lead free mark.
- The unleaded solder melts at a temperature about 40 °C higher than the ordinary solder, therefore, it is recommended to use the soldering iron having a temperature regulator.
- The ordinary soldering iron can be used but the iron tip has to be applied to the solder joint for a slightly longer time. The printed pattern (copper foil) may peel away if the heated tip is applied for too long, so be careful.

Section 2 Electrical Alignment

2-1. Procedure before and after Adjustment

Print out the record of user setting before performing the board replacement and adjustment. For the procedure of printing out the record, refer to Section 2-1-1.

If the history record of user setting cannot be printed due to the failure of this unit, reset the unit to the factory setting after performing the board replacement and adjustment. (Refer to Section 2-1-3.)

2-1-1. Procedure of Printing the User Setting Record (Before Adjustment)

Print out the top menu record of setting (last setting) and the list of all setting in the following procedure. After completing the board replacement and adjustment, if the setting is different from the contents of record printed before adjustment, adjust the setting so that it becomes the same as the setting before adjustment. Complete the procedure by adjusting the setting to the contents of top menu record.

Procedure

- Turn on the POWER.
 "READY" is displayed on LCD of this unit.
- 2. Press the OPEN button to open the door panel.
- 3. Set the thermal paper and close the door panel.
- 4. Press the jog dial.
- 5. Print out the record of the top menu (last setting) displayed on LCD.
- 6. Rotate the jog dial to display "MENU" on LCD, then press the jog dial.

All settings are printed on the thermal paper.

Note

Be sure to keep the printed sheets (user setting) because they are used as the reference to check the setting after adjustment.

2-1-2. User Setting Check Procedure (After Adjustment)

After completing the adjustment, reset this unit to the user setting in the following procedure.

Procedure

- 1. Print out the setting of this unit.
 (Refer to steps 1 to 6 of Section 2-1-1.)
- 2. Compare the setting list printed in step 1 with the setting list (used as reference) printed in step 6 of Section 2-1-1.

Note

If there is any difference between the settings, adjust to the reference setting using the jog dial. If there is no difference, press the jog dial. (Setting is completed.)

- After adjusting the setting to the reference setting, press the jog dial.
 - The adjusted setting is displayed on LCD.
- 4. Check that the contents of setting are correct, then press the jog dial.

2-1-3. Factory Setting

Procedure

1. Turn on the POWER while pressing the OPEN button.

Note

The beep sound is made immediately after turning on the POWER and after that again.

2. Release the OPEN button, then press the COPY button.

Note

Press the COPY button within 3 seconds after the second beep sound in step 1 is made.

- Check that this unit becomes the following state.(Change to the factory setting)
 - (1) Backlight is off.
 - (2) Backlight lights in green.
 - (3) Backlight lights in red.
 - (4) Flashes two times automatically.
 - (5) Beep sound is made.

2-2. Calibration and Electrical Conductivity Check

After replacing the KY-572 or MA-132 board, perform the calibration of brightness and contrast volumes, and perform the electrical conductivity check of each button.

1. Turn on the power while pressing the OPEN button and FEED button simultaneously.

The service mode starts.

Note

Do not release the OPEN button and FEED button before the beep sound is made.

- Press the jog dial two times. "SVC" and then "PATAN" are displayed on LCD in order.
- 3. Rotate the jog dial to display "BR:CAL".
- 4. Press the jog dial.

The calibration mode of brightness starts.

- (1) Rotate the BRIGHT volume fully to the left in the state that "BR:LEF" is displayed on LCD, then press the jog dial.
- (2) Rotate the jog dial to display "BR:RIG" and rotate the BRIGHT volume fully to the right in the state that "BR:RIG" is displayed on LCD, then press the jog dial.
- (3) Rotate the jog dial to display "BR:CEN" and return the BRIGHT volume to the center in the state that "BR:CEN" is displayed on LCD, then press the jog dial.
- 5. Rotate the jog dial to display "BACK", then press the jog dial.
- 6. Rotate the jog dial to display "CO:CAL".
- 7. Press the jog dial.

The calibration mode of contrast starts.

- (1) Rotate the CONT volume fully to the left in the state that "CO:LEF" is displayed on LCD, then press the jog dial.
- (2) Rotate the jog dial to display "CO:RIG" and rotate the CONT volume fully to the right in the state that "CO:RIG" is displayed on LCD, then press the jog dial.
- (3) Rotate the jog dial to display "CO:CEN" and return the CONT volume to the center in the state that "CO:CEN" is displayed on LCD, then press the jog dial.

Check the electrical conductivity of the COPY button.

Press the COPY button.

The error beep sound is made.

Note

The electrical conductivity of other buttons can be checked by calibration.

9. Turn off the power.


2-3. Head Voltage Adjustment

After replacing the thermal head or the switching regulator, perform this adjustment.

Note

Before starting the replacement, print out the stairstep test pattern which is used as the reference of adjustment.

Before replacing the switching regulator, measure the voltage of TP302 on the MA-132 board and take a note of the value.

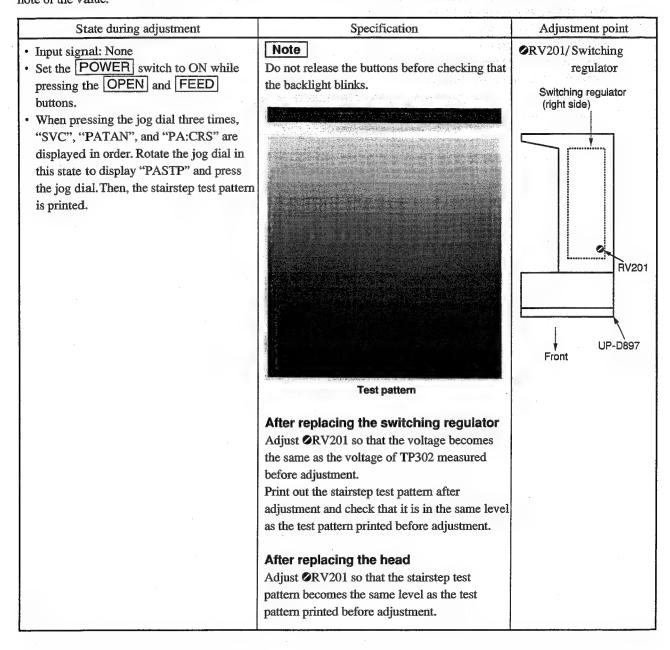

2-4. Initialization of Print Count History

After replacing the thermal head, perform the initialization of print count history.

State during adjustment	Specification	A	djustment point
Input signal: None		None	#1. %
Set the POWER switch to ON while			
pressing the OPEN and FEED buttons.			
When pressing the jog dial two times,			
"SVC" and "PATAN" are displayed in			
order. Rotate the jog dial in this state to			
display "RESET" and press the jog dial.	#*		
"R:MENU" is displayed and rotate the		. 777	
jog dial. After displaying the "R:H.PRN",	Check that "RESET" is displayed.		
then press the jog dial.	The initialization of print count history is completed (count becomes "0").		

2-5. Firmware Version Upgrade

After replacing the MA-132 board or IC103 on the MA-132 board, be sure to perform the firmware version upgrade.

Note

Never turn off the power during the firmware version upgrade.

Required equipment

 Personal computer (hereafter referred to as PC) (USB interface is mounted.)

OS: Windows 2000/XP

Driver: Driver software of UP-D897 is installed.

• Utility software for version upgrade

Note

Download the adjustment utility software from the URL described in the technical memo.

· Latest firmware

Note

Download the latest firmware from the GSP homepage or from the URL described in the technical memo.

USB cable

Note


Use the supplied USB cable or the equivalent (USB2.0 certified product).

Preparation

- 1. Connect this unit and PC using the USB cable.
- 2. Turn on the power of this unit and PC.
- 3. Install the driver software in PC. (Only when it is not installed.)
- 4. Copy the latest firmware to the directory same as the version upgrade utility.

Version upgrade

Start the version upgrade utility software.
 The version upgrade utility software screen is displayed.

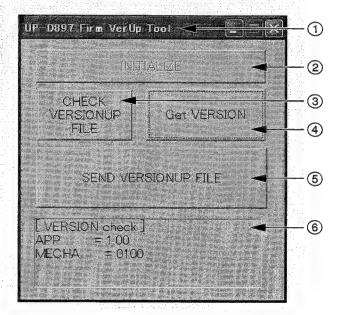

Number	Description
①	Model name to be upgraded
2	Search for the version upgrade model
3	Check of version upgrade file information
4	Check of firmware version written in this unit
(5)	Write the firmware in IC103 on the MA-132 board.
6	Status display window


2. Click the INITIALIZE button.

Note

In the case that the version upgrade utility software for the version upgrade of this unit is installed, the CHECK VERSIONUP FILE and Get VERSION buttons are enabled. In case that the cable is not connected or the version upgrade utility software for the version upgrade of this unit is not installed, the CHECK VERSIONUP FILE and Get VERSION buttons are disabled in gray out state.

When the CHECK VERSIONUP FILE and Get VERSION buttons are in gray out state, check the connection of PC and this unit, or the version upgrade utility software.

Click the Get VERSION button.
 The installed firmware version is displayed.

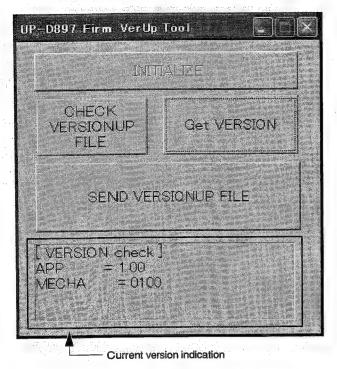

- 4. Check or take a note of the version displayed in step 3.
- Click the CHECK VERSIONUP FILE button.
 The VUP file information is displayed on another window.
- 6. Check the version and date to make sure that the downloaded firmware is the latest version.
- 7. Click the SEND VERSIONUP FILE button. The writing of firmware starts.

Note

When the writing starts, the LED backlight changes from staying lit to blinking, the beep sound is made. After this, the writing is completed when the backlight changes from blinking to staying lit.

- 8. Check that the writing is completed, and then turn the power of this unit off and on again (restart).
- 9. Start the version upgrade utility software.
- 10. Click the INITIALIZE button.
- 11. Click the Get VERSION button.
- 12. Make sure that the firmware is the latest version on the utility screen.
- 13. Exit the version upgrade utility software.

Procedure after completion of version upgrade

1. Turn off the power of this unit and PC, and then disconnect the USB cable.


Section 3 Circuit Description

3-1. System Control Block

3-1-1. Outline

System control is performed using IC104 (TMS320DM310ZHK22) on the MA-132 board. The system control block consists of the following blocks.

- CPU block (ARM925)
- · SDRAM control block
- External memory control block (Flash memory, etc.)
- DSP block (Image processing)
- External interface block (USB1.0, memory cards, video encoder, and CCD interface)

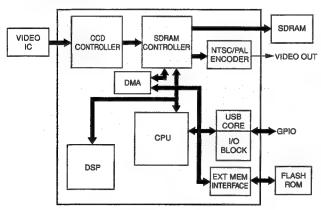

Block diagram of IC104

3-1-2. Block Diagram of Electrical Block

Refer to "Section 6 Block Diagram".

3-2. Head/Mechanical Control Block

3-2-1. Outline

The thermal head/mechanical control block is the circuit block of IC202 (CXD9182AGG) on the MA-132 board. There is an external I/O block except the thermal head/mechanical control block.

3-2-2. Platen Motor Control

A platen driving stepping motor controls the forward rotation/reverse rotation and rotation speed when IC202 controls driving transistors (Q603 to Q606).

IC202	Signal name	State			
Pin 46	PM_A1	L	L	Н	Н
Pin 47	PM_A2	Н	Н	L	L
Pin 92	PM_B1	L	Н	Н	L
Pin 93	PM_B2	Н	L	L	Н
			Forward r	otation	
			Reverse r	otation	

3-2-3. Thermal Head UP/DOWN Control

A thermal head UP/DOWN DC motor is driven using a driving circuit (IC201 on the SE-768 board). The DC motor can be rotated in the forward and reverse directions. It is controlled by IC202. Two types of optical position sensors (photo-interrupters PH201 and PH203 on the SE-768 board) detect the DC motor position and IC202 reads it.

Operation of Head UP/DOWN Motor

DC motor	IC202 (pin 1) or CN204 (pin 8) DCM_0	IC202 (pin 2) or CN204 (pin 7) DCM_1	Operation
Forward rotation	L	Н	The thermal head is raised.
Reverse rotation	Н	L	The thermal head is lowered.
Brake	Н	Н	Stop
Stop	L	L	Stop

State of position sensor

Position	IC202 (pin 53) or CN204 (pin 3) HEAD_P_SENSE	IC202 (pin 96) or CN204 (pin 5) HEAD_P2_SENSE	State
PRINTING	L	Н	Printing
HOME	Н	L	Standby (Usually)
DOOR OPEN	Н	Н	Door unlocked

3-2-4. Monitoring of Door Sensor

The door position is read using an optical door position sensor (photo-interrupter PH202 on the SE-768 board) by IC202.

State of door position sensor

Position	IC202 (pin 87) or CN204 (pin 2) DOOR_SENSE	State
CLOSE	Н	The door is closed.
OPEN	L 1	The door is opened.

3-2-5. Monitoring of Paper Sensor

Whether thermosensible paper is correctly set in this unit is detected using two pairs of optical paper sensors (phototransistors Q102 and Q302 on the KY-572 and SE-769 boards) and read using IC202.

State of paper sensor

IC202 (pin 43) or CN201 (pin 7) AD3 (PAPER_SENSE)	IC202 (pin 175) or CN205 (pin 2) AD5 (PP_EMP_SENSE)	State
Н	' H '	Paper
L	L	No paper
Н	L	No paper (Paper exists in an eject port and does not exist in a tray.)
L	Н	No paper (Paper does not exist in an eject port and exists in a tray.)

3-2-6. Monitoring of Head Temperature Sensor

The change in the resistance value of a thermistor in a thermal head is converted into a voltage and read using IC202. The A/D-converted voltage value corrects the density (gamma), controls the head cooling fan motor, and discriminates whether to clean the head.

3-2-7. Control of Head Fan Motor (for Head Cooling)

A head fan motor operates when IC202 controls driving transistors (Q607 and Q608). The head fan motor is turned on when the head temperature is more than approximately 62°C (cooling) or during printing.


Operation of head fan motor

IC202 (pin 94) HEAD_FAN	Operation
L STATE OF THE STA	OFF
H	ON

3-2-8. Read of Buttons

The OPEN, COPY, FEED, and jog dial buttons on the front panel are monitored. Each button (not including jog dial buttons UP and DOWN) is shifted to operation after the falling of a signal is detected. Jog dial buttons UP and DOWN are shifted to operation after the rising of a signal is detected.

IC202	Signal name	Function
Pin 7	KEY1	OPEN
Pin 98	KEY2	PRINT
Pin 55	KEY3	COPY
Pin 133	KEY4	FEED
Pin 6	JOG3	ENTER
Pin 54	JOG_UP	Rotates the jog dial upward.
Pin 97	JOG_DOWN	Rotates the jog dial downward.

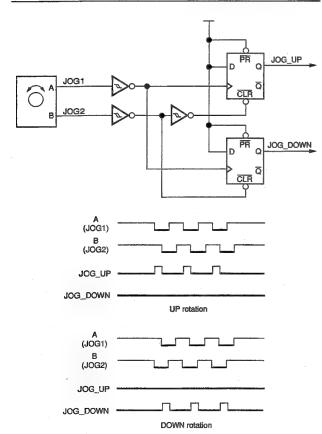


3-3. USB Interface

A USB2.0 high-speed interface is controlled by IC401 (UPD720122F) on the MA-132 board.


3-4. Thermal Head Block

3-4-1. Structure

A thermal head consists of 1280 dots (64 bits \times 20) per line. There are twenty pairs of blocks below.

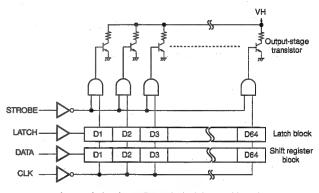

Internal circuit configuration of thermal head (corresponding to one port)

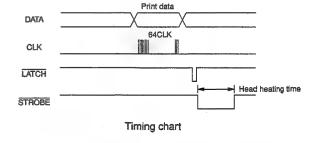

3-4-2. Basic Operation

Each signal is input from IC202 to the thermal head for operations below.

- (1) Print data is input to the shift register block in synchronization with a CLK pulse. (64-bit data)
- (2) The data input in step (1) is moved from the shift register block to the latch block when a latch pulse is input.
- (3) When a STB pulse is input, the "H" and "L" data of a latch block turn on and off an output-stage transistor and a resistor is heated for the color development of thermosensible paper.

3-4-3. Temperature Correction


The print energy required for thermosensible paper changes moment by moment due to the heating and thermal storage of a thermal head during change in room temperature or continuous printing. Therefore, the corresponding correction is required. In this unit, IC202 measures the temperature change of a thermal head from the internal thermistor of the thermal head. IC202 then converts the temperature change into 8-bit head temperature data and corrects the change in concentration for the temperature on which gamma characteristics were reflected. The concentration change can be corrected by controlling the width of a STB pulse.

3-4-4. Correction of Resistor Count


The total current flowing through a thermal head differs in the state where the resistors of the thermal head are all turn on or partially turned on. Therefore, an error occurs in the energy applied to each resistor. When print is made without correcting an error, a stripe occurs in the print result at the point where the number of resistors in which a thermal head is turned on changes rapidly. A circuit that corrects this stripe is incorporated into IC202.

Section 4 Troubleshooting


4-1. Print result is not satisfactory

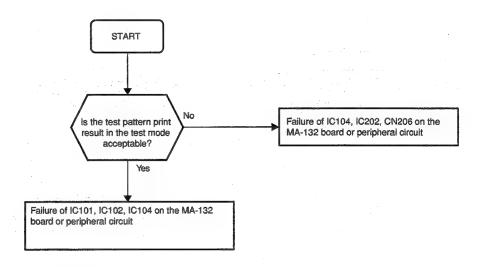


4-2. Print result density is too high or low

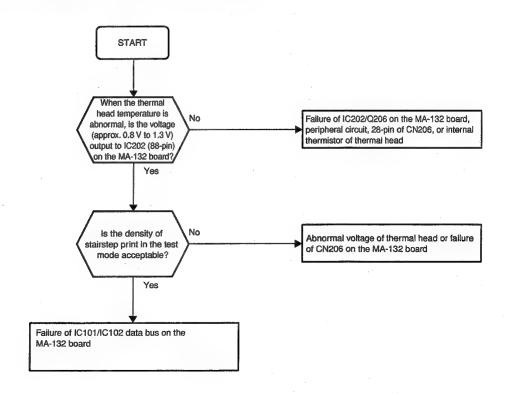


4-3. Trouble of determining presence or absence of paper

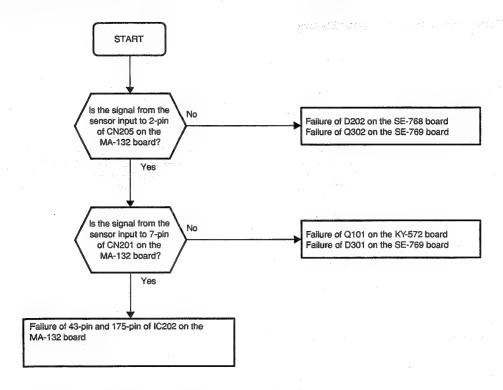


4-4. Thermal head UP/DOWN operation trouble

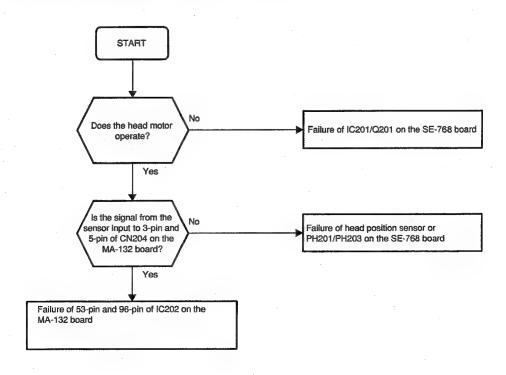

4-5. Feed operation trouble

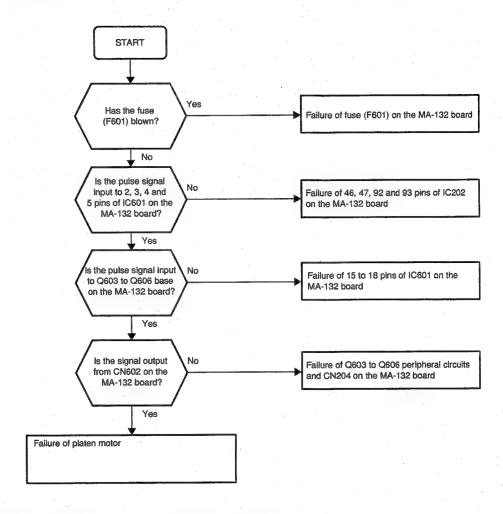

4-6. Trouble of determining door open/close

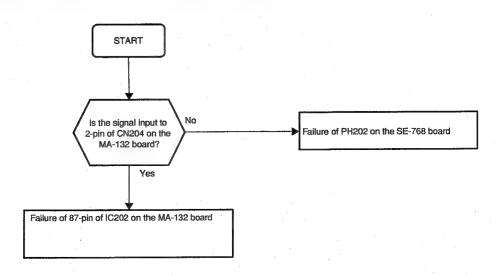

Section 5 Service Mode (Self-diagnosis Function)


5-1. Startup Procedure

- 1. Turn on the power while pressing the FEED button and OPEN button.
- Check that the backlight starts blinking after approximately 3 seconds, and then release each button. The service mode starts.

5-2. Service Mode Menu

When the unit is activated in service mode, the "SVC" menu is added one line above the "BEEP" menu. Press the jog dial in this "SVC" menu and rotate it up and down to display the item to be diagnosed, and then press it again. For the details of each item, refer to the following sections.

When the unit is activated in service mode, the time out function does not work. Therefore the unit does not exit the menu mode for approximately 20 seconds without any button operation.

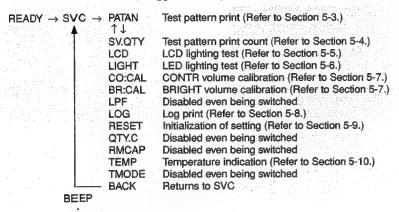

The usual menu continues.

5-3. Test Pattern Printing

Select the test pattern and press the jog dial to start the printing.


The preset number of test patterns are printed by setting the print count previously in "SV.QTY".


To interrupt the printing, press the FEED button or OPEN button.

When adjusting the head voltage, perform in the following conditions.


· Set "GAMMA" to "GA-2".


PATAN	PA:FF	All white	
	· •	Gray	
	PA:00	All black	
	PA:W/B	White in the left half, black	in the right half
	PA:B/W	Black in the left half, white	in the right half
	PA:FFh	All white (large size)	
	:	Gray (large size)	
	PA:00h	All black (large size)	
	PA:STP	Stairstep	
	PA:CRS	Cross step	
	BACK	Returns to PATAN	

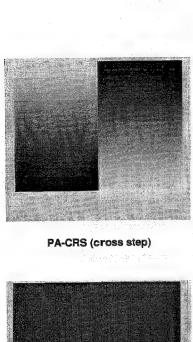


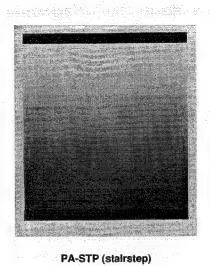




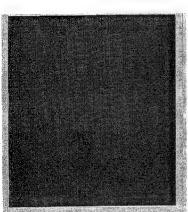

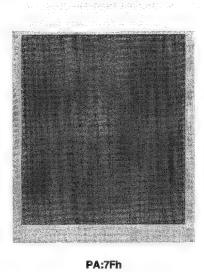


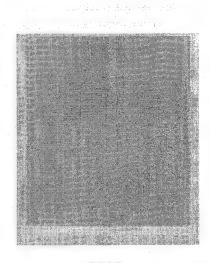





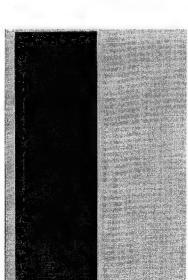



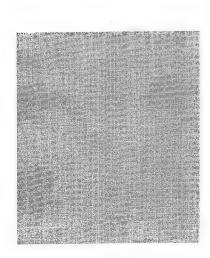

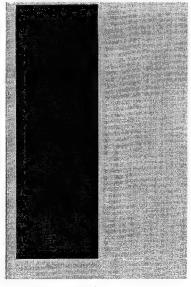


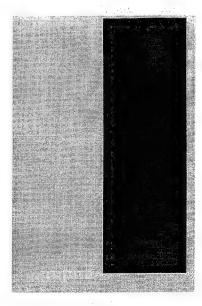


PA:FFh

PA:B/W

PA:W/B

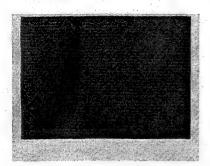


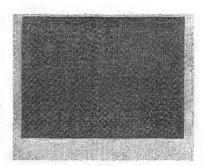


PA:00

PA:3F

PA:7F

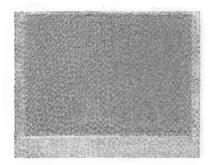

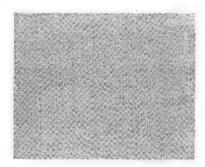

PA:BF

PA:FF

5-4. Test Pattern Print Count

Set the print count of the test pattern printing (PATAN).

SV.QTY SQ C

When setting the count to the continuous printing.

S.Q 9 When setting the count to "9".

S.Q:1

When setting the count to "1".

BACK Returns to SV.QTY.

5-5. LCD Lighting Check

Check the portion that is not lit by fully lighting the LCD.

LCD (full lighting) Check if it is lit.

Returns to LCD when the jog dial is

pressed.

BACK

Returns to LCD.



LCD (full lighting)

5-6. LED Lighting Check

Check if the backlight of the LCD and RECEIVING LED are lit.

LIGHT LI:RCV Check if the RECEIVING LED lights in green. (LCD lights in green.)

LI:AMB Check if the backlight of the LCD lights in amber.

LI:GRN Check if the backlight of the LCD lights in green.

BACK Returns to LIGHT.

5-7. Front Panel Volume Calibration

Perform the calibration of the CONTR volume and BRIGHT volume. Perform this calibration when replacing the flexible flat cable connecting the KY-572 board or MA-132 board. Each volume is adjustable at three points: left, center and right.

CO:CAL CO:RIG Press the jog dial in the state that the CONTR volume is fully rotated to the right.

CO:CEN Press the jog dial in the state that the CONTR volume is at the center click position.

CO:LEF Press the jog dial in the state that the CONTR volume is fully rotated to the left.

BACK Returns to CO:CAL.

BR:CAL BR:RIG Press the jog dial in the state that the BRIGHT volume is fully rotated to the right.

BR:CEN Press the jog dial in the state that the BRIGHT volume is at the center click position.

BR:LEF Press the jog dial in the state that the BRIGHT volume is fully rotated to the left.

BACK Returns to BR:CAL.

5-8. Log Printing

The print count or error log is printed. Set the printing paper in the paper tray.

LOG LO:OK The log is printed by pressing the jog dial.

BACK Returns to LOG.

Ver.1.00	lul 27 2005 - 0 - 1	Vilje N
T.PRN H.PRN T.ON F.ON H.TEMP 1.1501 2.1302 3.0000 4.0000 5.0000 6.0000 7.0000 8.0000	200 ② 100 ③ 100 ④ 10 ④ 10 ⑤ 30.7 ⑥ 50 ② 29.4 20 ② 27.9 0 0.0 0 0.0 0 0.0 0 0.0 0 0.0 0 0.0 0 0.0 0 0.0 0 0.0 0 0.0	

- ① Firmware version, date and time it is created
- ② T.PRN: Total print count
- (4) T.ON: Head total print count
- (4) T.ON: Total power ON time (unit: hour)
 (5) F.ON: Total fan ON time (unit: hour)
- 6 H.TEMP: Maximum head temperature (unit: °C)
 7 Error log: The eight error logs can be saved.
 - 1: The latest log, 8: The oldest log
 From the left: Error ID, total print count at
 the time of error, head temperature at the

time of error

Error ID

ID	Factor	
1102h	Head cannot be moved from home position to printing position. (During printing)	
1103h	Head cannot be moved from home position to printing position. (During feeding)	
1105h	Head cannot be moved from home position to printing position. (During cleaning)	
1204h	Head cannot be moved from home position to door open position. (During door open)	
1301h	Head cannot be moved from printing position to home position. (During mechanical initialization)	
1302h	Head cannot be moved from printing position to home position. (During printing)	
1303h	Head cannot be moved from printing position to home position. (During feeding)	
1305h	Head cannot be moved from printing position to home position. (During cleaning)	
1401h	Head cannot be moved from door open position to home position. (During mechanical initialization)	
1404h	Head cannot be moved from door open position to home position. (During door open)	
1501h	Head cannot be moved from arbitrary position to any of the positions (printing/home/door open). (During mechanical initialization)	
2100h	Thermistor is shorted. (Abnormally high temperature)	
2200h	Head is not connected. (Abnormally low temperature)	
2300h	Preheat time out	
2400h	Cool down time out	
3100h	Front paper sensor detects "no paper" during printing.	
5100h	Print pulse time out	
6101h	Prestart time out	
6102h	Paper feed before printing time out	

5-9. Initialization of Setting

RESET	R:ALL	Resets the unit to the factory
		setting.
	R:FAN	Resets the fan ON time.
		This is performed when replacing
		the fan.
	R.H.PRN	Resets the print count record to "0"
		This is performed when the head is
		replaced.
	R:MENU	Initializes the contents that are set
		in the menu.
	BACK	Returns to RESET.

5-10. Temperature Indication

The current temperature data is displayed. Example)

TEMP H.T:30.5 Head temperature: 30.5 °C BACK Returns to TEMP.

5-11. FEED Operation


When the FEED button is pressed in the state that the paper tray is open, the stepping motor is activated and the platen rotates in the paper eject direction.

5-12. Menu Lock

This is the function to prohibit the change of the setting from the menu and disable the switching of the front volume. It is disabled only when the unit is started normally.

Lock/Release Procedure

Turn on the power while pressing the jog dial. After checking that a beep sound is made after approximately 3 seconds, release the button. When the lock is activated, "LOCK" is displayed on the LCD and when the lock is released, "UNLOK" is displayed on the LCD. If the jog dial or the front volume is operated when the lock is activated, a "beep" alarm sound is made and "LOCK" is displayed on the LCD.



Section 6 Spare Parts

6-1. Notes on Repair Parts

1. Safety Related Components Warning

WARNING

Components marked \triangle are critical to safe operation. Therefore, specified parts should be used in the case of replacement.

2. Standardization of Parts

Some repair parts supplied by Sony differ from those used for the unit. These are because of parts commonality and improvement.

Parts list has the present standardized repair parts.

3. Stock of Parts

Parts marked with "o" at SP (Supply Code) column of the spare parts list may not be stocked. Therefore, the delivery date will be delayed.

4. Harness

Harnesses with no part number are not registered as spare parts.

In need of repair, get components shown in the list and repair using them.


5. Symbol

Protective earth (ground)




To identify any terminal which is intended for connection to an external conductor for protection against electric shock in case of a fault, or the terminal of a protective earth (ground) electrode.

Cover and Chassis Block

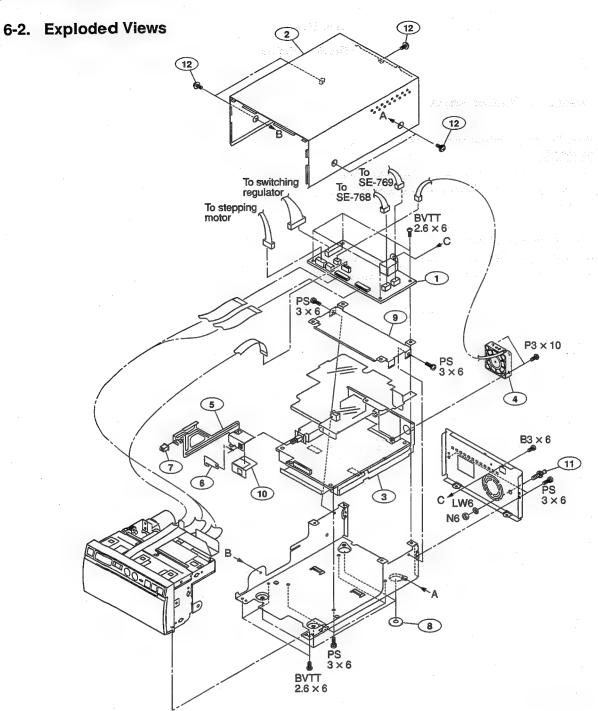

No.	Part No. SP Description	No.	Part No. SP Description
1 2 3	A-1078-885-A s MOUNTED CIRCUIT BOARD, MA-132(D) X-3698-532-3 s SUB ASSY, TOP COVER A 1-468-934-11 s REGULATOR, SWITCHING	12	4-886-821-11 s SCREW, M3X6 CASE (SILVER)
4	1-787-426-11 s FAN, DC (40 SQUARE)		7-623-425-07 s WASHER LW 6 (TYPE B) (for UC2/CED/SYN)
. 5	3-857-727-02 s ROD, POWER SW		7-682-149-04 s SCREW +P 3X10 (EP-FE/CU,NI,CR)
6	3-857-731-01 s STOPPER, SW ROD		7-682-547-04 s SCREW +B3X6 7-682-647-09 s SCREW +PS 3X6(EP-FE/ZNBK/CM2)
7	3-857-732-01 s BUTTON, POWER		7 -002-04/-U3 S DUKEW TPD DAO(EF-FE/ANDA/CHZ)
8	3-857-734-01 s FOOT		7-684-026-04 s NUT M6 TYPE2 (EP-FE/ZN/CM2)
9	3-857-735-02 s PLATE, SHIELD		(for UC2/CED/SYN)
10	3-857-736-01 s SHIELD, SW		7-685-862-09 s SCREW,+BVTT 2.6X6(EP-FE/ZNBK/CM2)
11	3-990-273-01 s TERMINAL, P.E. (for UC2/CED/SYN)		(-003-002-03 3 SCREW, TOVII 2.0AU (BE-FE) ANDRICHE)

Front Panel Bloc

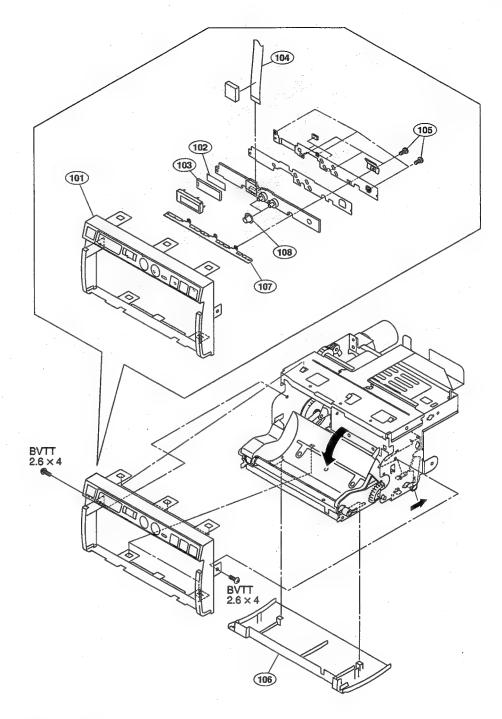


```
No. Part No. SP Description

101 X-3608-994-2 s FRONT PANEL ASSY(D)
102 1-780-223-11 s CONDUCTIVE BOARD, CONNECTION
103 1-805-786-11 s DISPLAY PANEL, LIQUID CRYSTAL
104 1-831-138-11 s CABLE, FLEXIBLE FLAT (24CORE)
105 3-713-791-45 s TAPPING SCREW M1.7

106 3-857-733-01 s PANEL, DOOR
107 3-863-111-02 s CUTTER
108 3-863-112-01 s KNOB, ROTARY
```


7-685-860-09 s SCREW, +BVTT 2.6X4 (EP-FE/ZNBK/CM2)

Mecha Deck Block 1

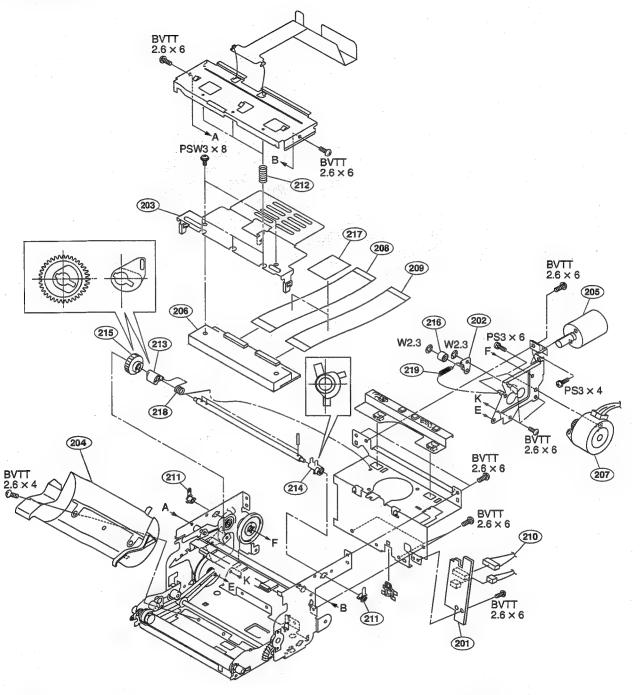


No.	Part No.	SP	Description	No.	Part	No.	SP	Description
201	A-1078-892-A	s	MOUNTED CIRCUIT BOARD, SE-768	213		7-753-0		
202	X-3608-983-1	S	SUB ASSY, TENSION ARM	214	3-85	7-754-0:	3 s	CÁM (R)
203	X-3608-987-3	S	ASSY, HEAT SINK	215	3-85	7-755-0	3 s	GEAR, CAM
204	X-3608-991-5	S	ASSY, PAPER HOLDER					
205	- X-3704-803-1	S	SUB ASSY, DC MOTOR (RP)	216	3-85	7-757-0:	l s	ROLLER, TENSION
200				217	3-85	7-783-0	4 s	SUPPORT (FFC)
206	1-479-084-11	S	HEAD, THERMAL (LVE6426SS)	218	3-98	6-503-03	s	SPRING, TORSION
207	1-787-289-11	S	MOTOR, STEPPING	219	3-98	6-932-0	l s	SPRING, MOTOR TENSION (1N) S
208	1-830-214-21	S	CABLE, FLEXIBLE FLAT (26 CORE)					
209			CABLE, FLEXIBLE FLAT (28 CORE)					
210	1-963-347-12	S	HARNESS, MA-SE768		3-66	9-596-01) s	WASHER, 2.3 (PLA)
			·		7-68	2-645-04	1 s	ISO SCREW+PS3X4 CR-N (13)
211	3-857-741-01	s	BEARING		7-68	2-647-09	s	SCREW +PS 3X6(EP-FE/ZNBK/CM2)
212	3-857-749 - 01	S	SPRING, HEAD (10N)		7-68	2-948-01	8	SCREW +PSW 3X8

Mecha Deck Block

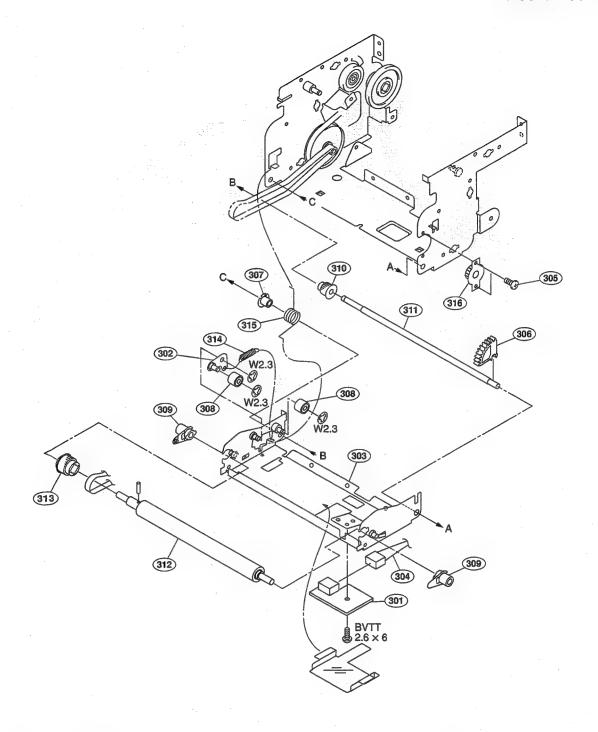

No.	Part No. SP	Description	No.	Part No.	SP Description
303 304	X-3608-983-1 s X-3608-984-1 s 1-963-348-12 s	MOUNTED CIRCUIT BOARD, SE-769 SUB ASSY, TENSION ARM SUB ASSY, MD DOOR HARNESS, MA-SE769 SCREW (M2X4), LOCK ACE, P2	313 314	3-857-762-01 3-857-763-03 3-857-770-01	s SHAFT, FULCRUM s PLATEN s PULLEY, PLATEN s SPRING, DOOR EXTENSION (4N) s SPRING, DOOR OPEN (C)
307 308 309	3-857-757-01 s	SPACER, FULCRUM ROLLER, TENSION BEARING, PLATEN		3-669-596-00	s DAMPER, OIL s WASHER, 2.3 (PLA) s SCREW, +BVTT 2.6X6 (EP-FE/ZNBK/CM2)

Mecha Deck Block 3

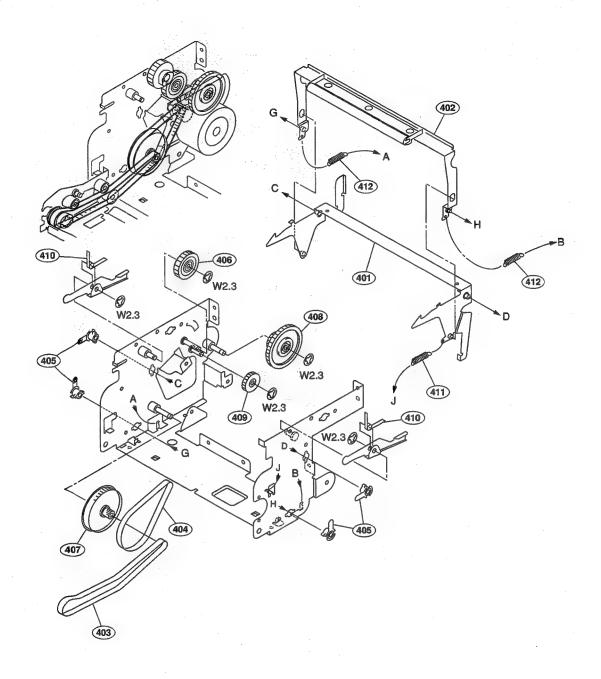

No.	Part No.	SP	Description	No.	Part No.	SP	Description
401 402 403 404 405	X-3608-986-2 3-854-457-01	. S . S	SUB ASSY, DOOR LOCK SUB ASSY, PINCH ARM TIMING BELT (220TN10-5.0T) TIMMING BELT (140TN10-4.0T) BEARING	411 412	3-990-805-03 3-669-596-00	. a	SPRING, LOCK (1.3N) S SPRING, PINCH (7.5N) WASHER, 2.3 (PLA) SCREW, +BVTT 2.6X4 (EP-FE/ZNBK/CM2)
406 407 408 409 410	3-857-750-03 3-857-751-03	S	PULLEY, IDLER				

6-3. Electrical Parts List

KY-572(D)G BOARD	MA-132 (D) G BOARD
Ref. No. or Q'ty Part No. SP Description	Ref. No. or Q'ty Part No. SP Description
C101 1-162-970-11 s CAPACITOR CERAMIC 0.01MF/25V B	1pc A-1159-429-A s MOUNTED CIRCUIT BOARD, MA-132(D)G
C102 1-162-970-11 s CAPACITOR CERAMIC 0.01MF/25V B C103 1-162-970-11 s CAPACITOR CERAMIC 0.01MF/25V B C105 1-162-970-11 s CAPACITOR CERAMIC 0.01MF/25V B	BZ201 1-544-886-11 s BUZZER, PIEZOELECTRIC
C105 1-162-970-11 s CAPACITOR CERAMIC 0.01MF/25V B C106 1-162-970-11 s CAPACITOR CERAMIC 0.01MF/25V B	C101 1-125-777-11 s CAPACITOR CERAMIC 0.1MF/10V C102 1-125-777-11 s CAPACITOR CERAMIC 0.1MF/10V
C107 1-162-970-11 s CAPACITOR CERAMIC 0.01MF/25V B C108 1-165-128-11 s CAPACITOR CERAMIC 0.22MF/16V F C109 1-165-128-11 s CAPACITOR CERAMIC 0.22MF/16V F C110 1-165-128-11 s CAPACITOR CERAMIC 0.22MF/16V F	C103 1-125-777-11 s CAPACITOR CERAMIC 0.1MF/10V C104 1-125-777-11 s CAPACITOR CERAMIC 0.1MF/10V C105 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF
CN101 1-817-368-61 s CONNECTOR, FFC/FPC (ZIF) 24P	C106 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C107 1-165-989-11 s CAPACITOR, CERAMIC 10MF (2012)
D101 6-501-137-01 s DIODE CL-375TD/SYG-D-TS D102 8-719-064-52 s DIODE CL-191YG-CD-T	C108 1-125-777-11 s CAPACITOR CERAMIC 0.1MF/10V C110 1-165-989-11 s CAPACITOR, CERAMIC 10MF (2012) C111 1-125-777-11 s CAPACITOR CERAMIC 0.1MF/10V
EN101 1-477-089-31 s ENCODER (ROTARY)	C112 1-165-989-11 s CAPACITOR, CERAMIC 10MF (2012) C113 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF
IC101 8-759-465-98.s IC BU9728AKV-E2	C114 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C115 1-125-777-11 s CAPACITOR CERAMIC 0.1MF/10V
Q102 6-550-988-01 s TRANSISTOR CPT-184S-C-TS-BCD	C116 1-125-777-11 s CAPACITOR CERAMIC 0.1MF/10V
R102 1-216-864-11 s CONDUCTOR, CHIP (1608) R105 1-216-853-11 s RESISTOR, CHIP 470K 1/16W (1608) R106 1-216-864-11 s CONDUCTOR, CHIP (1608) R107 1-216-864-11 s CONDUCTOR, CHIP (1608) R109 1-216-864-11 s CONDUCTOR, CHIP (1608)	C117 1-125-777-11 s CAPACITOR CERAMIC 0.1MF/10V C118 1-125-777-11 s CAPACITOR CERAMIC 0.1MF/10V C120 1-125-777-11 s CAPACITOR CERAMIC 0.1MF/10V C121 1-125-777-11 s CAPACITOR CERAMIC 0.1MF/10V C123 1-125-777-11 s CAPACITOR CERAMIC 0.1MF/10V
RV101 1-227-680-12 s RESISTOR, VAR, CARBON 10K RV102 1-227-680-12 s RESISTOR, VAR, CARBON 10K	C124 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C127 1-125-777-11 s CAPACITOR CERAMIC 0.1MF/10V C128 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF
\$101 1-771-884-31 s SWITCH, TACTILE \$102 1-771-884-31 s SWITCH, TACTILE \$103 1-771-884-31 s SWITCH, TACTILE	C129 1-125-837-91 s CAPACITOR, CHIP CERAMIC1MF/6.3V C130 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF
VDR101 1-801-924-21 s VARISTOR, CHIP (1608) VDR102 1-801-924-21 s VARISTOR, CHIP (1608) VDR103 1-801-924-21 s VARISTOR, CHIP (1608)	C132 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C133 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C134 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C135 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C136 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF
	C137 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C138 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C139 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C140 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C141 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF
	C142 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C143 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C144 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C145 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C147 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF
	C149 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C150 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C152 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C153 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C155 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF
	C156 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C157 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C158 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C159 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C160 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF
	C161 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C162 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C163 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C164 1-164-943-11 s CAPACITOR, CHIP CERAMIC 0.01MF C166 1-125-777-11 s CAPACITOR CERAMIC 0.1MF/10V

6-7

R204

R205

1-218-889-11 s RESISTOR, CHIP 56K 1/10W (1608)

R104

R105

1-208-643-11 s RESISTOR CHIP 22 1/16W (1005)

1-208-643-11 s RESISTOR CHIP 22 1/16W (1005)

6-11


6-4. Supplied Accessories

```
*1: [UP-D897(UC2)]
*2: [UP-D897(CED)]

Ref. No.
or Q'ty Part No. SP Description

1pc A *2 1-551-631-00 s POWER-SUPPLY CORD (IEC)
1pc A *1 1-556-813-22 s CORD, POWER
1pc 1-824-211-41 s CORD, CONNECTION
1pc 3-623-865-01 s SHEET, HEAD CLEANING (ACLYLIC)
1pc 3-863-306-03 s COMPACT DISC
```

Section 7
Block Diagram

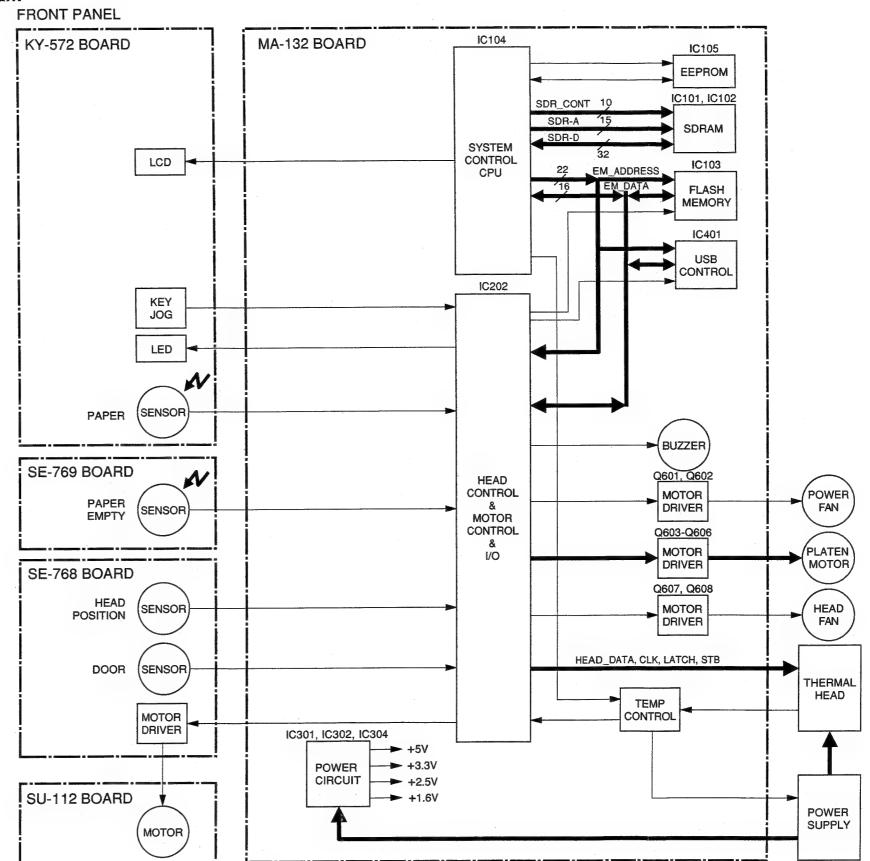


Section 8 Schematic Diagrams

Index

Board name	Page			
KY-572	8-2			
MA-132	8-3			
SE-768	8-9			
SE-769	8-9			
SU-112	8-9			
Frame Wiring	8-10			


8-1

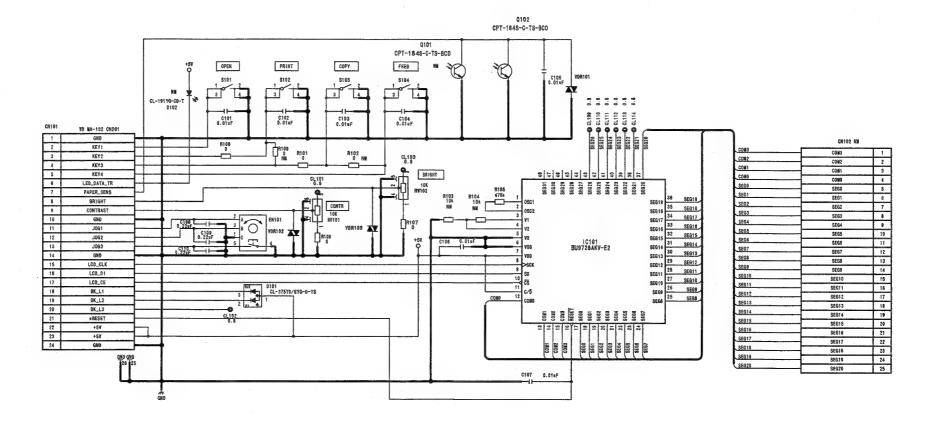


KY-572

D

C

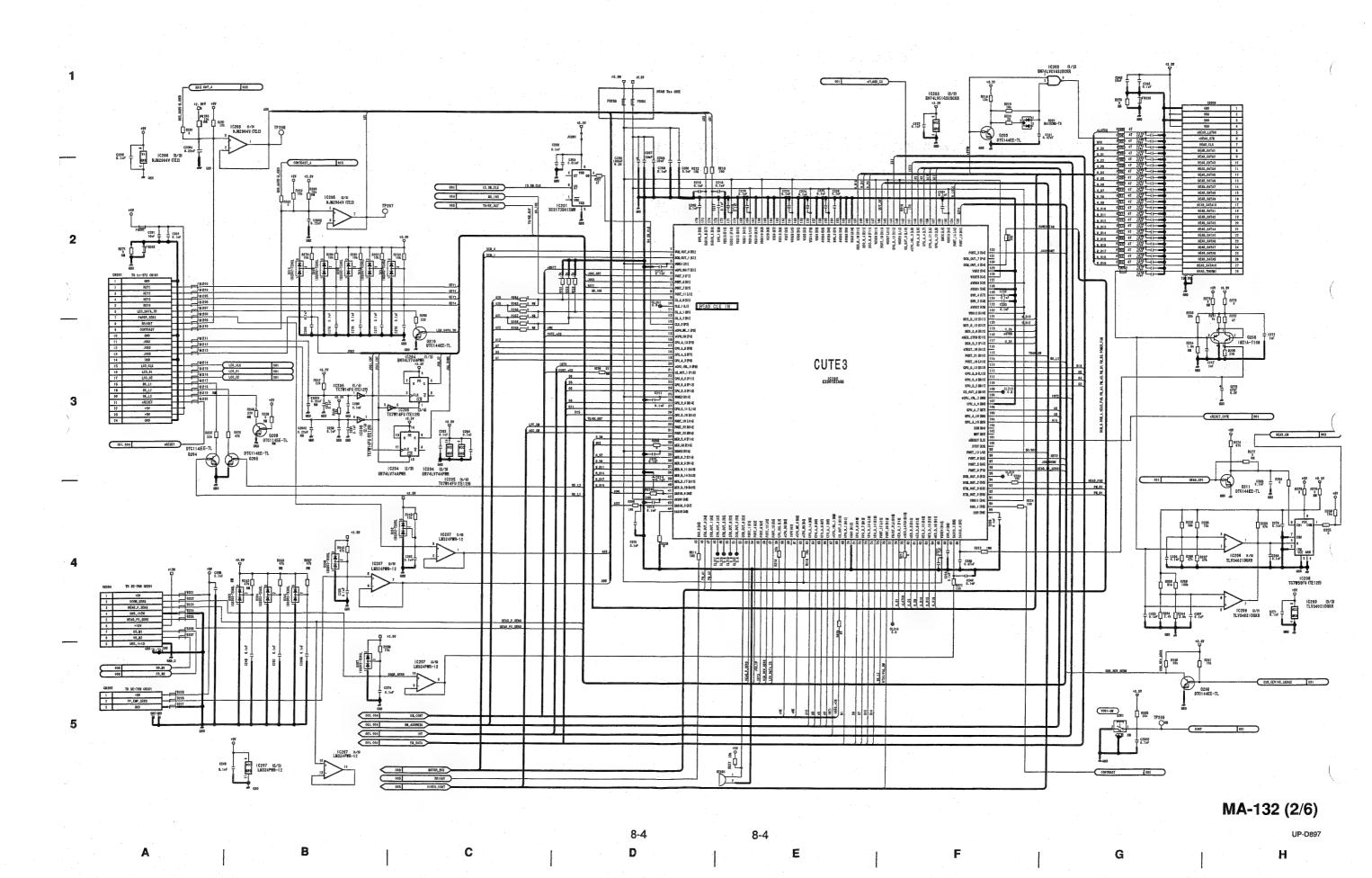

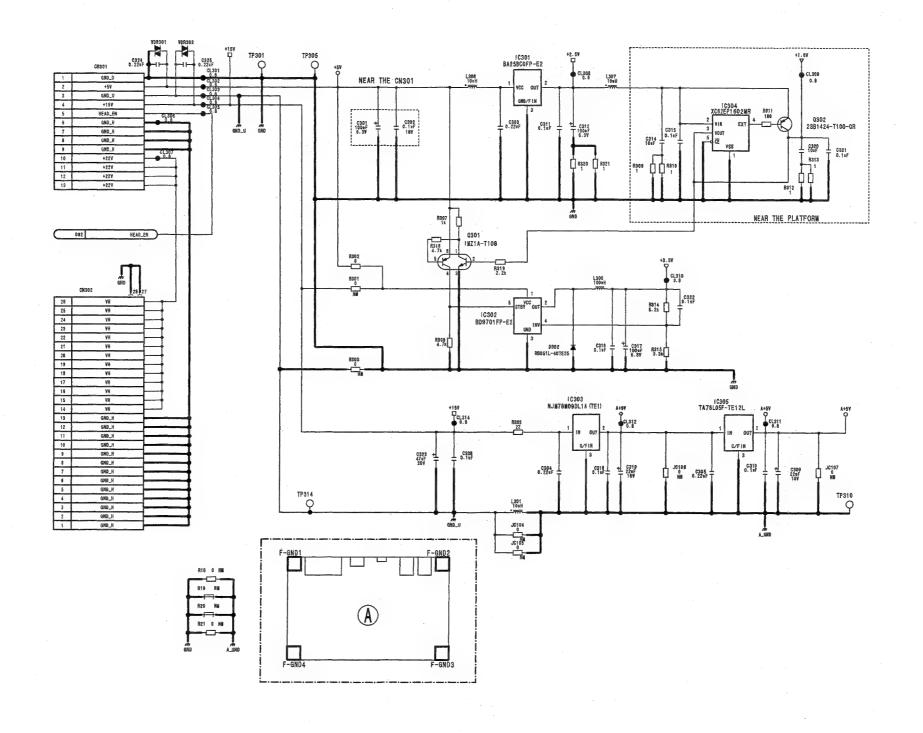

MA-132 (3/6)

1

С


8-5 D 8-5


F

G

G

Н

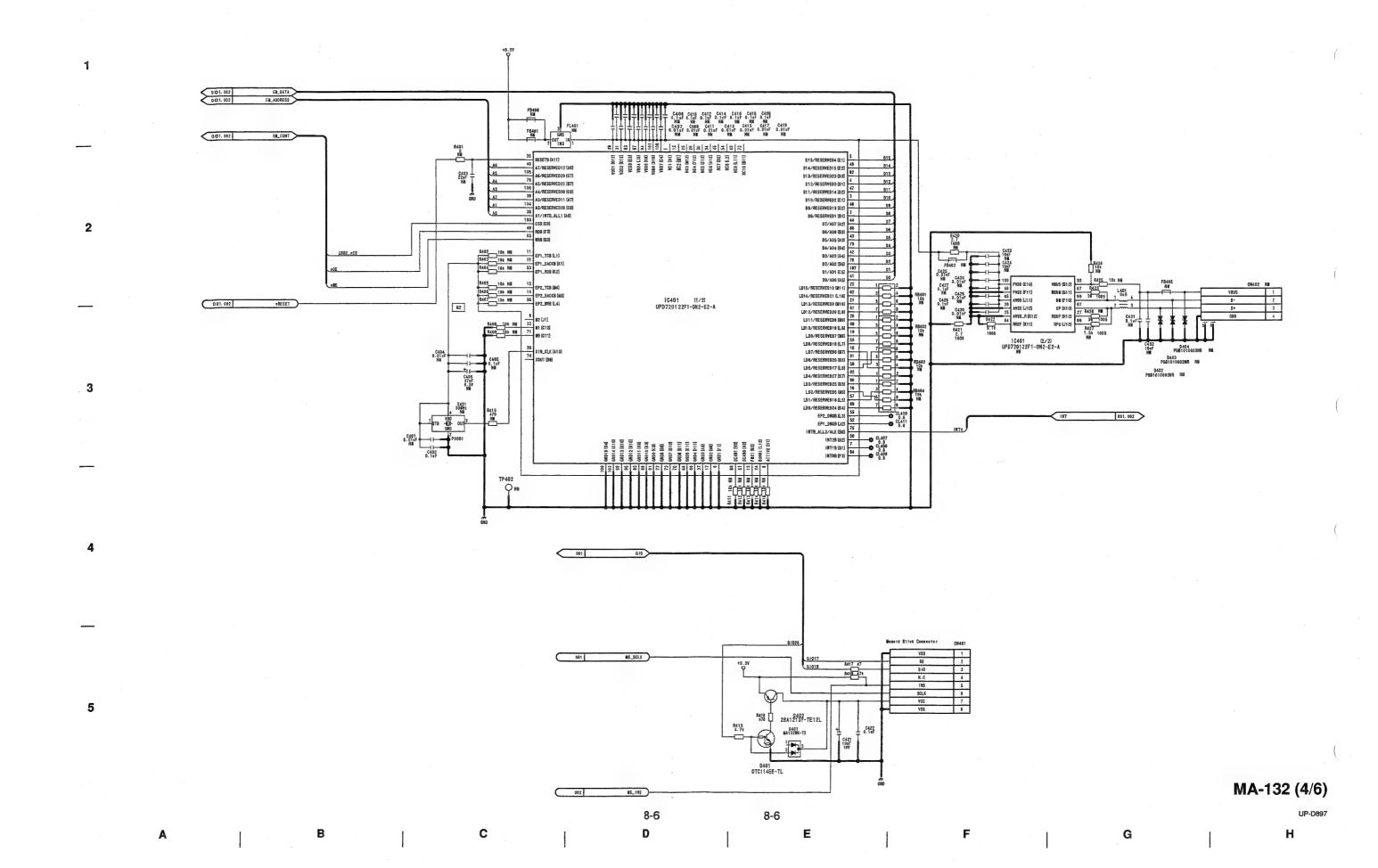

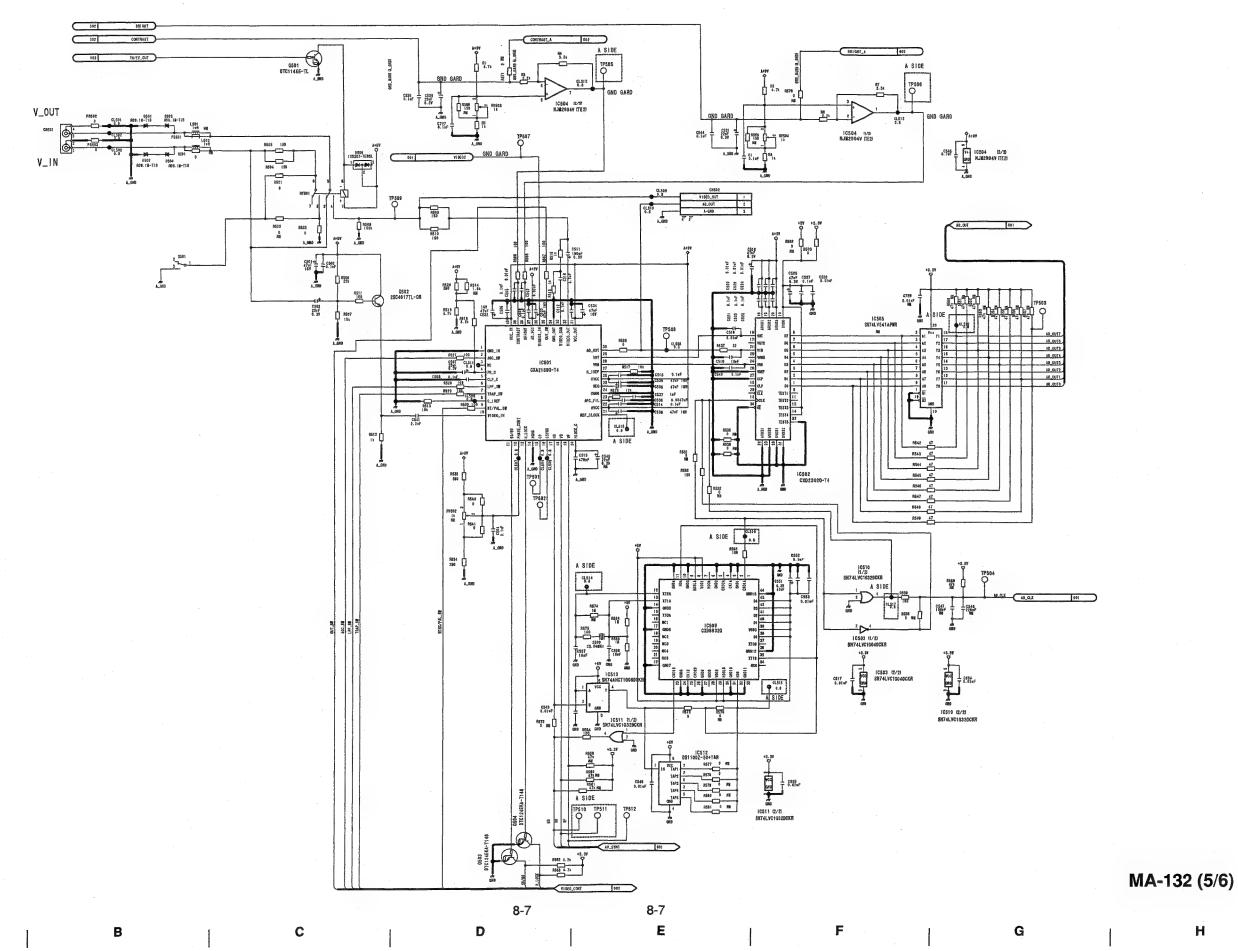

F400/125V/1.25A A609 560 R604 R625 2.2k R802 Q607 2SD999-T1-CLCK Q605 2SD992-Z-E2 PATTERN WIDETH =1.5mm Q608 280992-Z-E2 PATTERN WIDETH =1.5mm
PATTERN WIDETH =1.5mm
PATTERN WIDETH =1.5mm
PATTERN WIDETH =1.5mm
PATTERN WIDETH =1.5mm R601 | R627 PLACE OUT OF THE BLOCK R807 R610 100 R628 FIG12 R813 10k R814 R617 | R616 | R619 | R620 | R621 | R622 | 47k | 47k | 47k | PM_B2 PM_B2 PM_A2 PM_B1 POWER_FAN
HEAD_FAN
DCM_0
DCM_1

2


3

5


MA-132 (6/6)

UP-D897


8-8 8-8 B C D E F G H

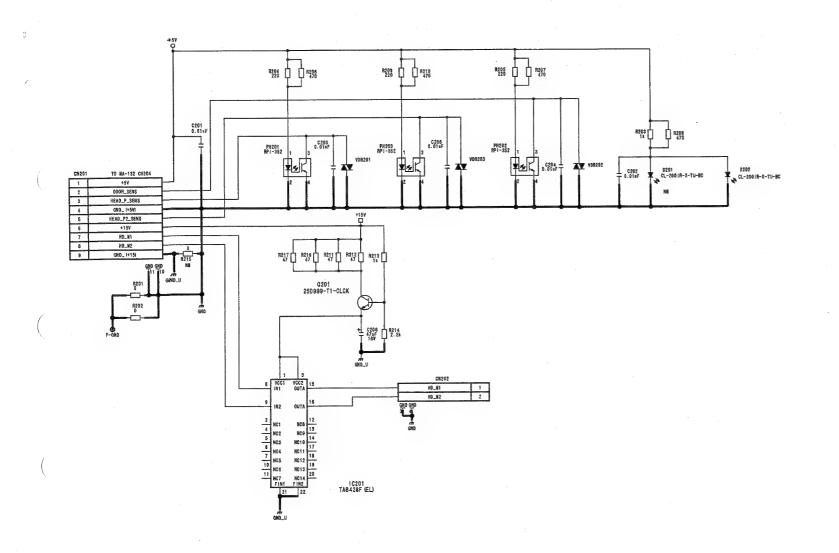

SE-768

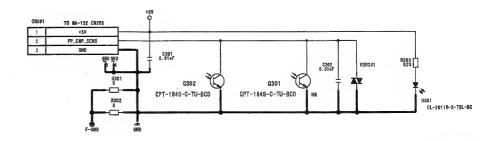

SE-769

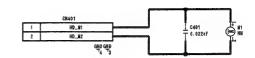


G

SU-112

8-9 8-9 1 **B** 1 **C** 1 **D** 1 **E**

UP-D897


н

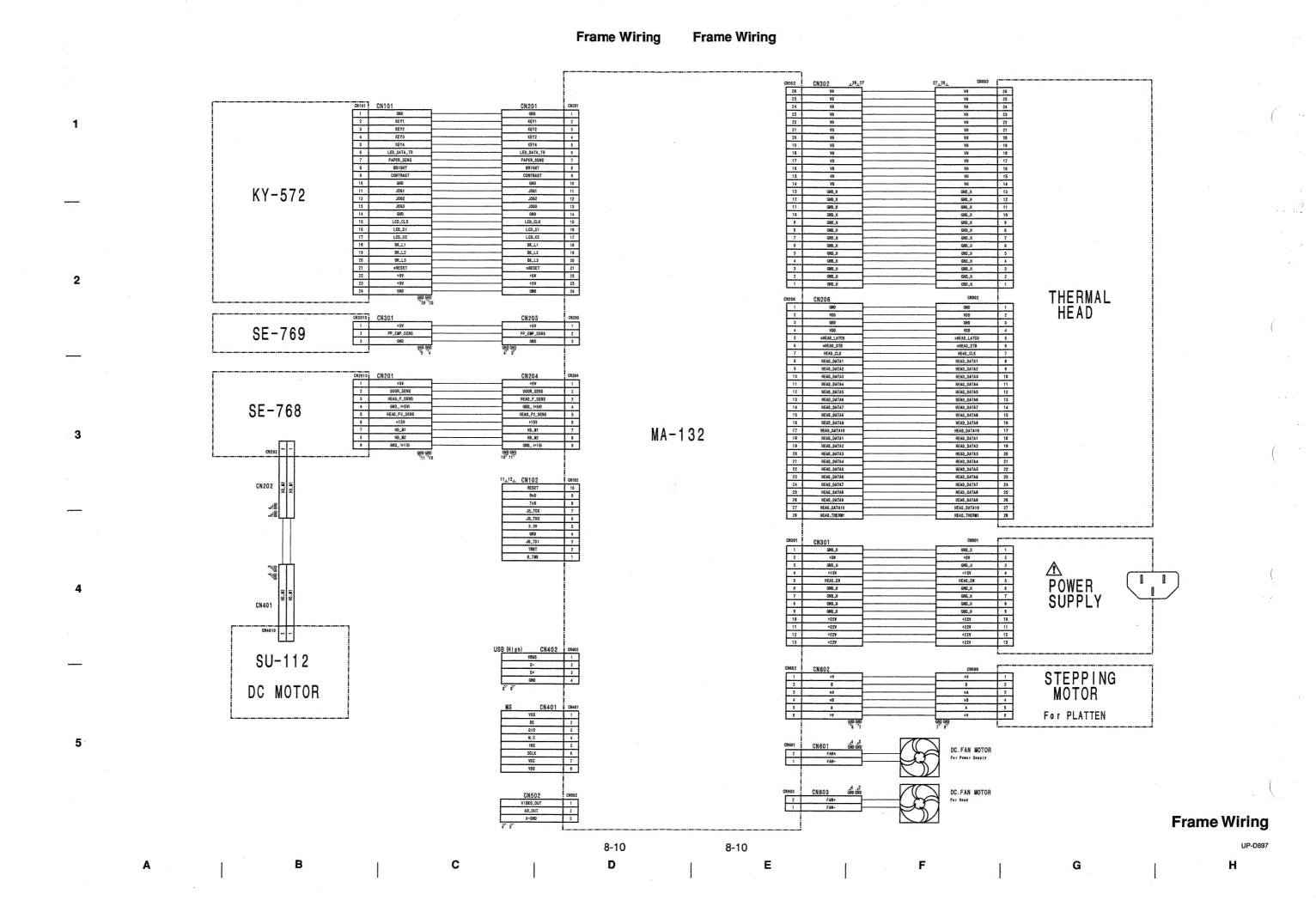


Section 9 Board Layouts

Index

Board name	Page			
KY-572	9-1			
SE-768	9-1			
SE-769	9-1			
SU-112	9-1			
MA-132	9-2			

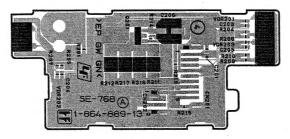

SE-768 -A SIDE-SUFFIX: -13

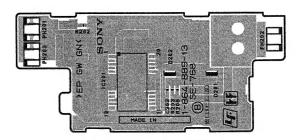

SE-768 -B SIDE-SUFFIX: -13



KY-572 -A SIDE-SUFFIX: -12



SE-769 -A SIDE-SUFFIX: -13




SE-769 -B SIDE-SUFFIX: -13




KY-572 -B SIDE-SUFFIX: -12




SU-112 -A SIDE-SUFFIX: -13

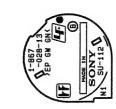

SU-112 -B SIDE-SUFFIX: -13

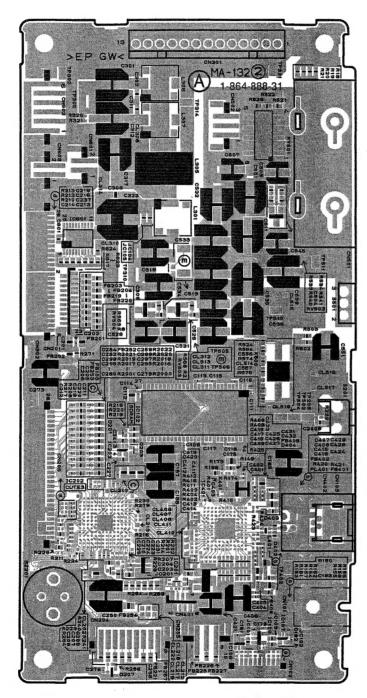

MA-132 -A SIDE-SUFFIX: -31

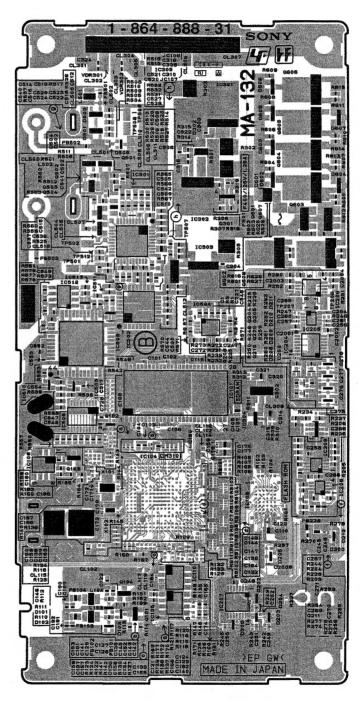

MA-132 -B SIDE-SUFFIX: -31

SAFETY CHECK-OUT

After correcting the original service problem, perform the following safety checks before releasing the set to the customer:

Check the metal trim, "metallized" knobs, screws, and all other exposed metal parts for AC leakage. Check leakage as described below.

LEAKAGE TEST

The AC leakage from any exposed metal part to earth ground and from all exposed metal parts to any exposed metal part having a return to chassis, must not exceed 0.5 mA. Leakage current can be measured by any one of three methods.

- A commercial leakage tester, such as the Simpson 229 or RCA WT-540A. Follow the manufacturers' instructions to use these instruments.
- A battery-operated AC milliammeter. The Data Precision 245 digital multimeter is suitable for this job.
- Measuring the voltage drop across a resistor by means of a VOM or battery-operated AC voltmeter. The "limit" indication is 0.5 V, so analog meters must have an accurate lowvoltage scale. The Simpson 250 and Sanwa SH-63Trd are examples of a passive VOM that is suitable. (See Fig. A)

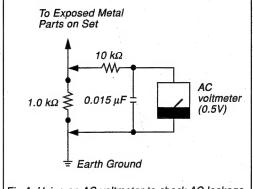


Fig A. Using an AC voltmeter to check AC leakage.